


How to pair your hearing aids to your iPhone


Questions? Call 847-662-9300

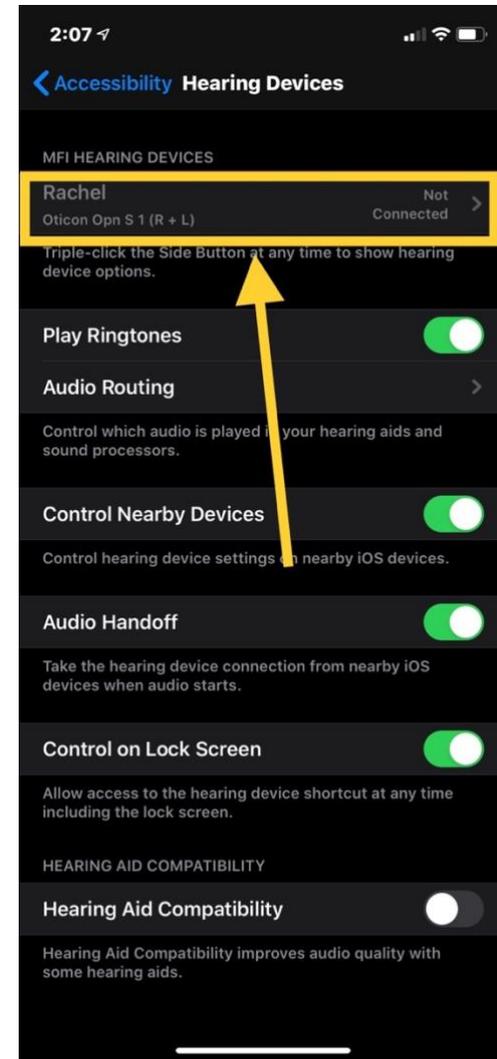


1 Open the Settings app on your phone. Scroll down to "Accessibility".

2 Scroll down to "Hearing Devices".


3 If your hearing aids were paired to your device before, you'll see this with your name. Click this. If it says "Searching...", skip to step 5.

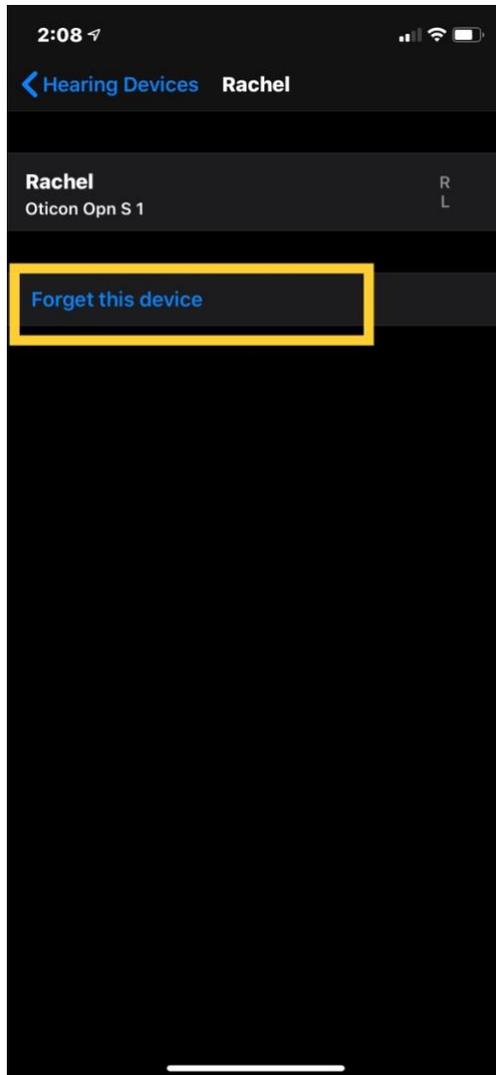

4 Click "Forget this Device" so that your hearing aids can re-pair.

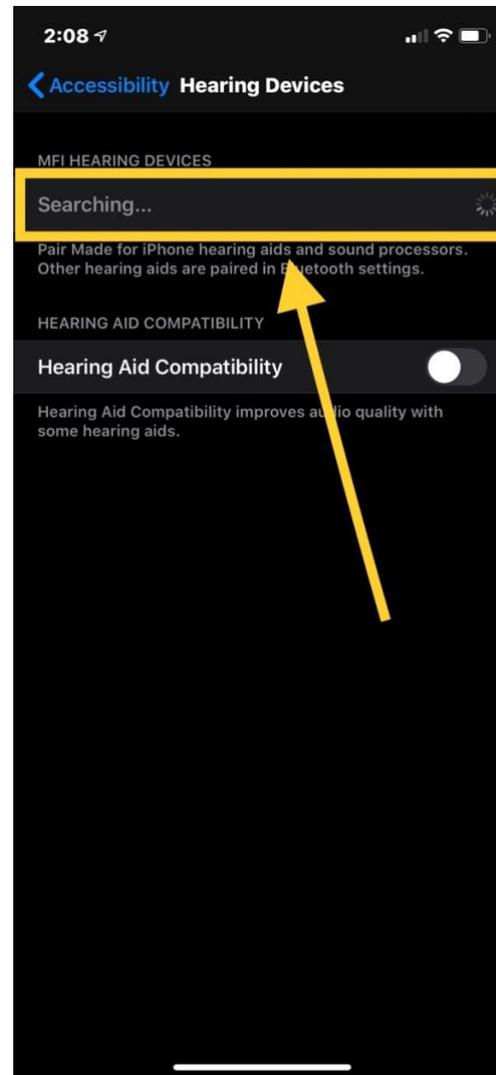

5 Now you'll see that your phone is searching for your hearing aids.

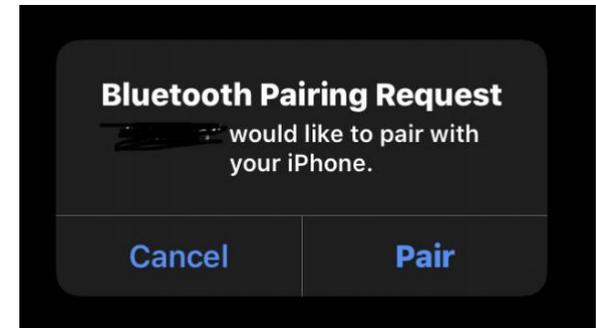

6 For disposable batteries: Open the battery door(s) and then close them again.
For rechargeable batteries: Either put in charger and take out, or turn them off by holding the down button for 4 seconds until it turns red, and then holding the down button for 6 seconds until it turns green again.

By restarting your hearing aids, they will be in pairing mode for 2 minutes. Your hearing aids should appear on the "searching" screen during this time.

Click the hearing aids when they appear and then click "Pair to Hearing Aids" (shown above). Note that if you have 2 hearing aids, you'll have to click "Pair" two times.

You are now paired and ready to use your hearing aids with your phone! Enjoy!